

T-Top 10 Autoclave

- ◆ Low Cost - Top Performance
- ◆ Class S Sterilizer (Closed-Door Drying)
- ◆ T-Connect: Advanced Remote Diagnostics
- ◆ Large Chamber Capacity
- ◆ Short Cycles



P/N: T-Top 10S 110V



Plug and Play

No dedicated installation connections required

- ◆ Multicolour user-friendly touch screen with intuitive icons
- ◆ Easy - just start and run programs
- ◆ Compact and ergonomic, designed to fit into your space



Reliable Performance

Quality product that will do all you need and expect

- ◆ Easy to operate, maintain and use
- ◆ Internal memory - storing up to 500 cycles
- ◆ Advanced technology to deliver



Traceability & Connectivity

Advanced Technology that provides real-time data intelligence

- ◆ 2 USB Ports
- ◆ Wi-Fi connection & T-Connect software
- ◆ Remotely access cycle history, get notifications and alerts.
- ◆ Paperless documentation



Innovation

Best in class technology

- ◆ Programming and customization functions
- ◆ Specially designed user interface to speed up and optimize daily workload
- ◆ Time saver setting to set tests and sterilization cycles to start at any time








Cloud based paperless software

T-Top devices feature the exclusive and technologically advanced T-Connect, a cloud based IoT solution for your sterilization center.




T-Top Sterilization Cycle Programs

Icon	Program name	Sterilization Times (w/ Drying)*	Load Type	Type of Use
	Unwrapped 273°F	25 min	Unwrapped Instruments	For immediate use
	Wrapped 273°F	56 min	Wrapped Pouches	For storage
	Unwrapped 250°F	46 min	Unwrapped Delicate Instruments	For immediate use
	Wrapped 250°F	72 min	Wrapped Delicate Instruments	For storage
	Wrapped 250°F / 2 Kg	42 min	Small Loads of up to 4.4 lbs.	For quick turnaround of small loads

* Cycle times above are an average for a standard load and could vary based on load weight, type, and external conditions

T-Top Technical Specifications

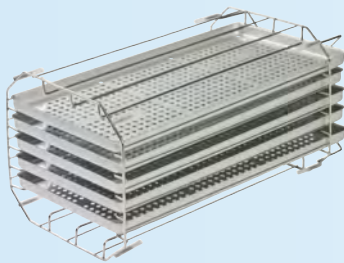

Chamber size	10 (W) x 18 (D)
Volume	21L
Overall dimensions	23" (D) x 18.1" (W) x 18" (H)
Max # of trays	5 (sold with 4 trays, additional available for purchased)
Tray dimensions	14.7" (D) x 7.3" (W) x 0.6" (H)
Weight	99 lbs.
Water filling	manual

Tuttnauer USA Co.

345 Oser Ave.
Hauppauge, NY 11788
Tel: +800 624 5836, +631 737 4850
Fax: +631 737 0720
E-mail: info@tuttnauerUSA.com

Visit Our Website:

www.tuttnauer.com/us

